Take a Close Look at Skin, Hair & Scalp

images with polarized light

Visioscope® PC 35









hair facial erythrosis

- USB camera with polarisation filters for images with & without reflection of the skin surface
- showing skin structure as well as pigmentation & blood flow
- freezing the image
- visualization and evaluation in the CSI-software
 - desquamation (dry skin) with Corneofix®
 - sebum production with Sebufix®
 - pigmentation spots
 - skin aging/ sun damage
 - lines and wrinkles
 - pores
 - skin impurities (reddening, acne, comedones, etc.)
 - hair structure, dandruff and condition of scalp
 - and many more applications




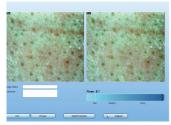

evaluation pores

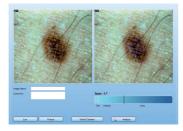

evaluation spots

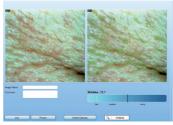

evaluation wrinkles



Technical Data:

Camera works with Windows® 7/8/10; Interface/Power supply: USB 2.0; Dimensions: approx. 12 x 5.5 x 5.5 cm; Cable length 1.5 m; Illumination: 2 x 8 white light LEDs; Area of view: 10 x 8 mm; Resolution: 1280 x 1024; Filters: 2 polarisation filters, one switchable; Camera button to switch the polarisation filter and freeze the image

